

A photograph showing a medical professional using a Henso Wound Closure Device on a patient's arm. The device consists of a yellow mesh strip with white adhesive tabs and a central red line. A pair of surgical forceps is shown pulling one of the white tabs, which is attached to the mesh. The background is a light blue gradient.

Wound Closure Device

Product Manuals

Henso Medical (Hangzhou) Co., Ltd.


Non-invasive suture on the skin surface, so that the wound healing beautiful and no trace



About Henso Hold

HENSO Hold Wound Closure Device is indicated to close acute wounds such as lacerations and surgical incisions. It's a new non-invasive technique that can reduce tissue trauma, improve patient comfort and possibly cause less scarring after healing.

The device consists of medical adhesive tape, a pair of polypropylene hasps and a liner. The tape is a porous nonwoven or PU film polypropylene hasp that is designed to allow for a non-invasive closure of a wound and to encourage a tight closure of the wound for an early and successful healing. Wound Closure Device is individually packed and terminally.

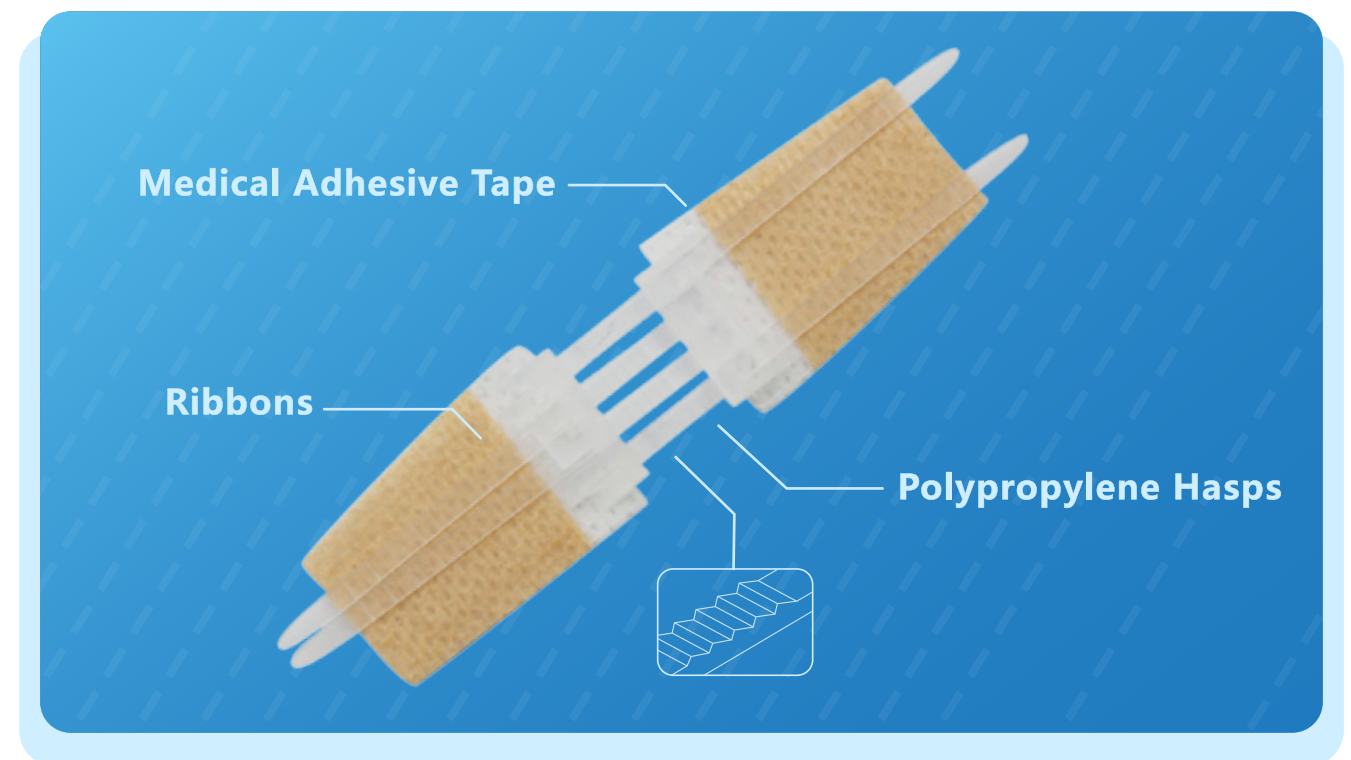

Ribbons


- Precise and unique structure design ensures the convenience operation and sufficient closing force when pulling.

Polypropylene Hasps

- Self-locked structure design ensures the ribbon will not be loosened during the treatment.
- It can be tightened gradually, especially for the edema wounds.
- Made of special poly material, provides strong closure force.

Medical Adhesive Tape

- Air-permeable and perfect biocompatibility ensures the device placed on human skin for long period.
- Strong and Stable holding power provides enough strength to close the wound.
- 180° peel strength ensures peeling easily without any glue residue.

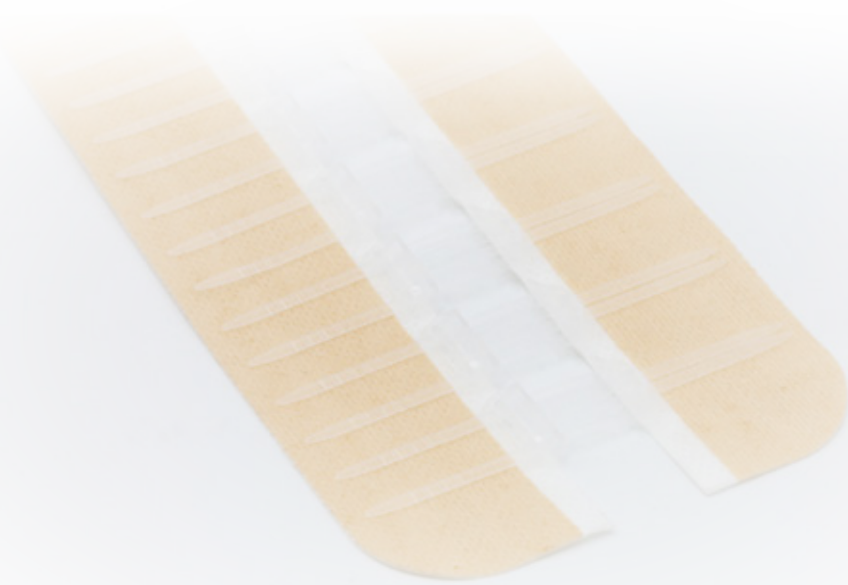

Applicable Cases

- Obstetrics and Gynecology Surgery
- Gastrointestinal Surgery
- Cardiothoracic Surgery
- Hepatobiliary Surgery
- Orthopedics Surgery
- Maxillofacial Surgery
- Emergency Surgery
- Pediatric Surgery
- Normal Wounds




Advantage

- No suture, less scars
- Less pain for the patients
- Quick to apply, easy to use
- Easy to clean the wound skin
- Ideal for Accident & Emergency
- Reduce the wound infection rate
- Adjustable and good for wound closure


Operation Instruction

1. Stop bleeding and clean the wound.
2. Open the package. ①
3. Take the wound closure device off the paper. ②
4. Stick on the wound skin. ③
5. Press the adhesive ribbons carefully. ⑤
6. Apply enough devices according to the wound degree (each device should be 2mm away from others). ④
7. Tighten the hasps gently. ⑤
8. Cut the long ribbons, apply dressing if possible. ⑥⑦
9. Remove the devices lightly after recovery. ⑧

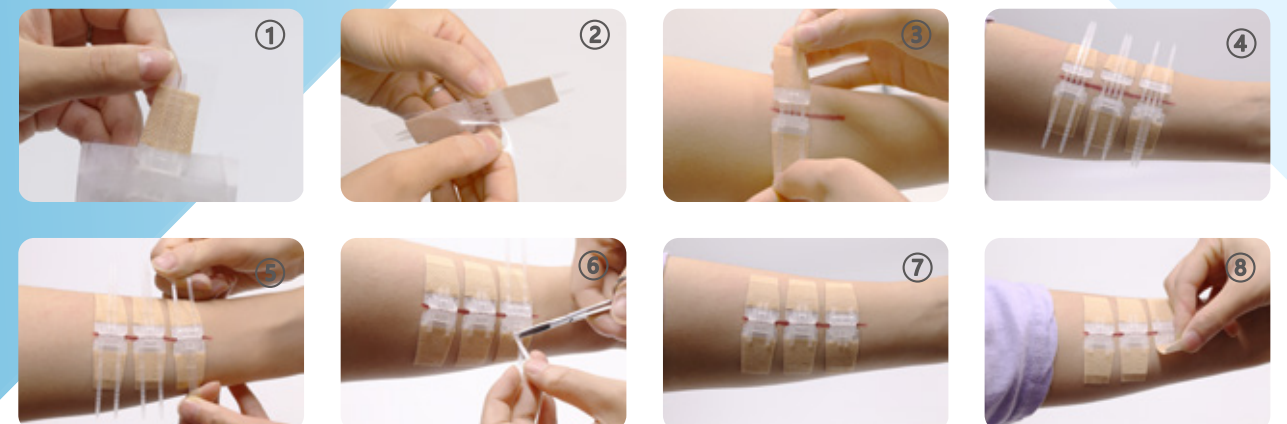



Single

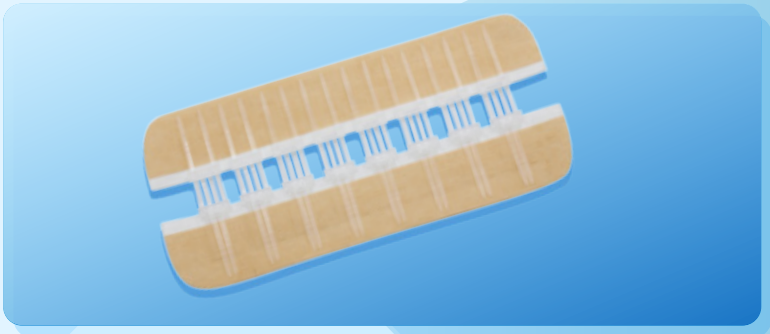

Combination

Specification

Item	Specs	Pack	Materials	Color
Single	PWB 1001	1pc/sheet,10sheet/box	Non-woven	White
	PWB 1002	2pc/sheet,10sheet/box		
	PWB 1003	3pc/sheet,10sheet/box		
	PWF 1001	1pc/sheet,10sheet/box		Skin
	PWF 1002	2pc/sheet,10sheet/box		
	PWF 1003	3pc/sheet,10sheet/box		
	PWM 1001	1pc/sheet,10sheet/box	PU film	
	PWM 1002	2pc/sheet,10sheet/box		
	PWM 1003	3pc/sheet,10sheet/box		

Item	Specs	Pack	Materials	Color
Combination	PZB 1003	10cm/sheet,10sheet/box	Non-woven	White
	PZB 1004	13cm/sheet,10sheet/box		
	PZB 1005	15cm/sheet,10sheet/box		
	PZB 1006	18cm/sheet,10sheet/box		
	PZB 1007	20cm/sheet,10sheet/box		
	PZB 1008	23cm/sheet 10sheet/box		
	PZB 1009	25cm/sheet,10sheet/box		
	PZB 10010	28cm/sheet,10sheet/box		
	PZB 1012	33cm/sheet,10sheet/box		
	PZF 1003	10cm/sheet,10sheet/box		
	PZF 1004	13cm/sheet,10sheet/box		
	PZF 1005	15cm/sheet,10sheet/box		
	PZF 1006	18cm/sheet,10sheet/box		
	PZF 1007	20cm/sheet,10sheet/box		
	PZF 1008	23cm/sheet,10sheet/box		
	PZF 1009	25cm/sheet,10sheet/box		
	PZF 1010	28cm/sheet,10sheet/box		
	PZF 1012	33cm/sheet,10sheet/box		
	PZM 1003	10cm/sheet,10sheet/box	PU film	Skin
	PZM 1004	13cm/sheet,10sheet/box		
	PZM 1005	15cm/sheet,10sheet/box		
	PZM 1006	18cm/sheet,10sheet/box		
	PZM 1007	20cm/sheet,10sheet/box		
	PZM 1008	23cm/sheet,10sheet/box		
	PZM 1009	25cm/sheet,10sheet/box		
	PZM 1010	28cm/sheet,10sheet/box		
	PZM 1012	33cm/sheet,10sheet/box		